



F-900D WHITE-SNOW PRINCESS VAPORIZER (MACHANISM)

1. Function:

Vapor and Sterilize skin

The vapor steam consists of rich and various active hydronium ion, can stimulate blood circulation and supply the water element in the cells, nourish and moisten the skin, coinstantaneous can sterilizer skin, when massage to use it can get perfect effect.

2. Structure:

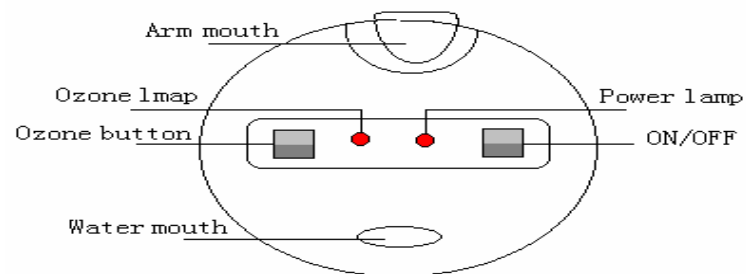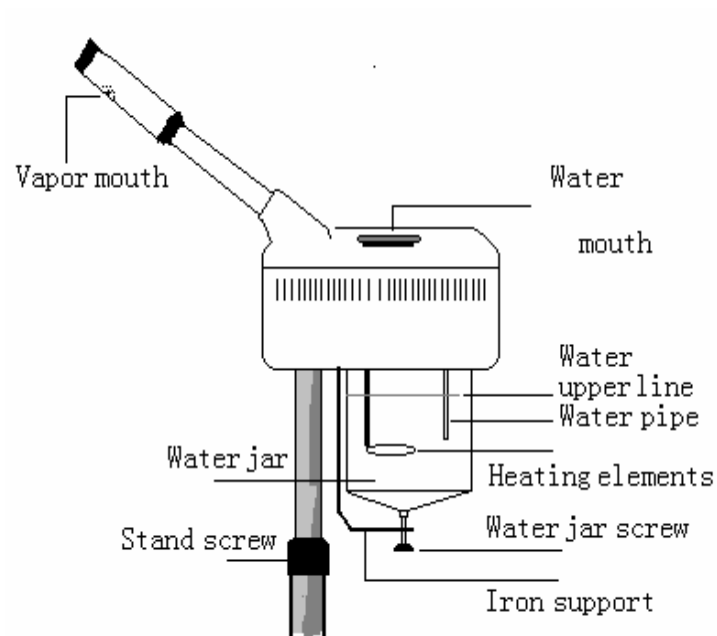

3. Instruction:

1. Dip distilled water into water glass through water mouth but don't over the water upper mark line always.
2. Press ON/OFF then the power indicator lights, then the machine is ready for work.
3. Press ozone button then the ozone indicator lights, ozone function is ready for work.

4. When the water in the water glass is under the water low mark line, and the machine will cut off the power automatic, after water 5 minutes later, the machine can continue work.
5. Machine should be put a suitable place where the vapor steam can freely blow from client's chin to each part of face.
6. The machine should keep earth well.

4. Caution:

1. Only when the machine can spray out vapor steam well, the operation begins.
2. Don't over 25 minutes every operation, after over the operation, clean water glass always.
3. If found the deposit on the heating element, dip the heating element in the vinegar acid 5 minutes and can remove it to avoid the deposit growing and prolong the heating element use.

5. Operating Time Table

Skin conditions	Space	Operating time
Oily, Poor blood circulation	20 cm	15 minutes
Normal	25 cm	10 minutes
Sensitive, Dry	30 cm	5 minutes
Telangiectasis, Ozone-violet Dermatitis	30 cm	5 minutes

6. Specification:

Model No.	Volt.	Power	Hz.	Temp.	Size	G.W.
F-900D	110/220V	750W	50/60Hz.	5-40°C	93*33*28cm	8.5kg