

F-908 HERB VAPORIZER MANUAL

1. Function:



1. Vapor and Sterilize skin

The vapor steamer consists of rich and various active hydronium ion, can stimulate blood circulation and supply the water element in the cells, nourish and moisten the skin, coinstantaneous can sterilizer skin, when massage to use it can get perfect effect.

2. Herb vapor

Use proper herb, have a great effect to treat comedo, telangiectasis, dermatitis, whiten and nourish skin.

2. Structure:

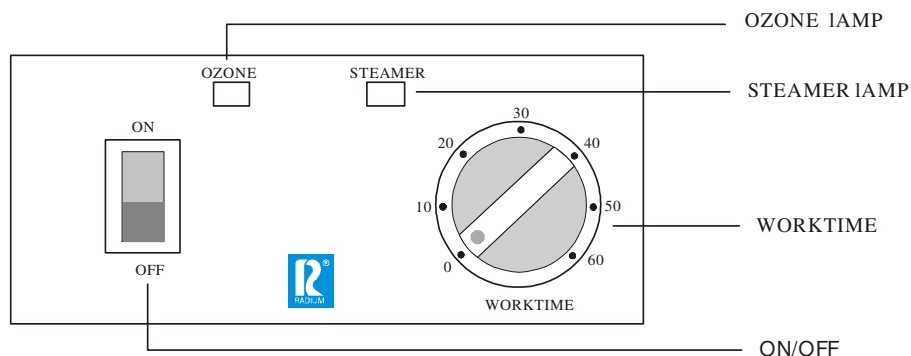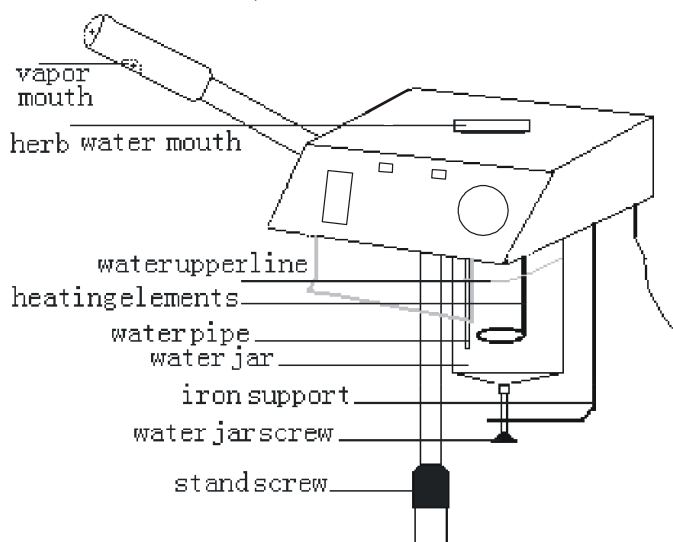

3. Instruction:

1. Open the herb cover and dip the distilled water into water glass through water mouth but don't over the water upper mark line always. At this moment, the machine can

work to operate normal vapor steam. (note: should tight the water glass by adjust screw to avoid the vapor steam leaving out from the water glass.)

2. Operate herb vapor, open the herb cover and put proper herb into the filter cup and tight the herb cover.
3. Turn on the timer then the machine is ready for work.
4. Press ozone button then the ozone indicator lights, ozone function is ready for work.
5. When the water in the water glass is under the water low mark line, and the machine will cut off the power automatic, after water 5 minutes later, the machine can continue work.
6. Machine should be put a suitable place where the vapor steamer can freely blow from client's chin to each part of face.
7. When the timer is over, the machine stops spray vapor and returns to wait state.
8. The machine should keep earth well.

4. Caution:

1. Only when the machine can spray out vapor steamer well, the operation begins.
2. Don't over 25 minutes use herb vapor function every operation, after over the operation, clean water glass always.
3. When herb vapor function works, don't coinstantaneous use Ozone function to avoid infect the herb vapor effect.
4. In the operation, when open the herb cover, first press the button (START/PAUSE) to avoid scald by vapor steam.
5. If found the deposit on the heating element, dip the heating element in the vinegar acid 5 minutes and can remove it to avoid the deposit growing and prolong the heating element use.

5. Operating Time Table

Skin conditions	Space	Operating time
Oily, Poor blood circulation	20 cm	15 minutes
Normal	25 cm	10 minutes
Sensitive, Dry	30 cm	5 minutes
Telangiectasis, Ozone-violet Dermatitis	30 cm	5 minutes

6. Specification:

Model No.	Volt.	Power	Hz.	Temp.	Size	G.W.
F-908	110/220V	750W	50/60Hz.	5-40°C	90*36*36.5cm	11.0kg